




The fundamental basis of Nuclear Magnetic Resonance (NMR) measurements on fluid-bearing rocks is that the decay or relaxation time of the NMR signals (T_2) is directly related to the pore size. The NMR signal detected from a fluid-bearing rock therefore contains T_2 components from every different pore size in the measured volume. Using a mathematical process known as inversion, these components can be extracted from the total NMR signal to form a T_2 spectrum or T_2 distribution, which is effectively a pore size distribution. From this distribution, various petrophysical parameters such as porosity, permeability, and free and bound fluid ratios can be measured or inferred.

However, as might be expected, things are not quite so simple in practice. The T_2 distributions calculated from NMR signals depend on many other parameters besides simply the sizes of the fluid-containing pores. Other factors that influence the shape of the distributions include: the nature of the fluid (in particular the viscosity); the spacing between the data points (echoes) in the NMR signal; the logging speed; the wettability of the rock; and the hydrogen index of the fluid(s) being observed. The variability of these factors in the natural environment can lead to errors in interpreting the logging information collected during the exploration process.

In the laboratory, most of these additional factors can be eliminated or at least controlled by making measurements on fully saturated core samples and on ones that are at irreducible water saturation. These measurements help to distinguish between fluid and pore size effects. In the logs, these effects have a significant influence on the recorded data making interpretation difficult. Calibration of logs refers to comparing data from the logs with data from well characterised core plugs taken from the same well, in an attempt to make interpretation of the log data easier and more accurate. The value of this process is recognised by logging experts:

“...careful calibration of petrophysical interpretation models with core data is essential for obtaining accurate permeability and bound fluid answers from NMR logs, as compared to using default parameters in the petrophysical model.”

Schlumberger leaflet on NMR core analysis

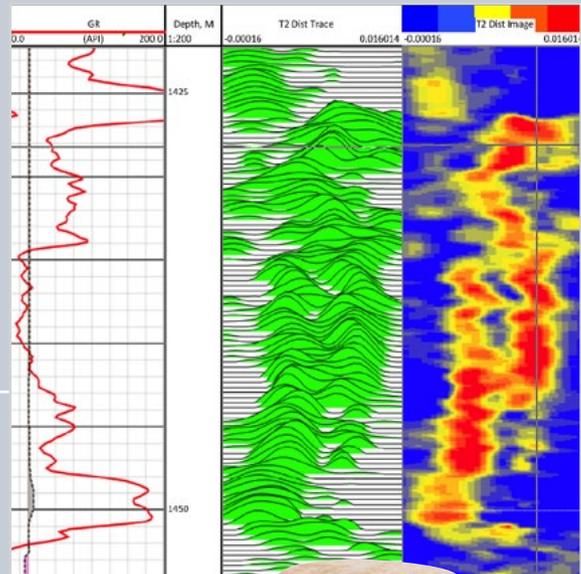

The Business of Science®

NMR log calibration from laboratory core measurements



What factors can be calibrated?

The two NMR log measurements most often calibrated using core plugs are the T_2 cut-off dividing the bound and free fluid, and permeability modelling coefficients.

T_2 cut-off

The so-called “ T_2 cut-off” in a T_2 distribution is the T_2 value that divides the small pores that are unlikely to be producible from the larger pores that are likely to contain free fluid. The integral of the distribution above the T_2 cut-off is a measure of the free fluid (mobile fluid) in the rock, and is clearly influenced by the position of the T_2 cut-off point, as shown in Figure 1. The portion of the curve below the cut-off is known as bound fluid and is made up of the clay bound fluid and the capillary bound fluid.

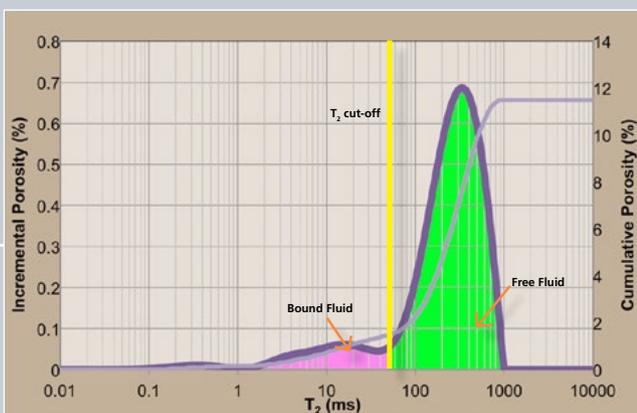


Figure 1 - T_2 distribution showing integral for free fluid



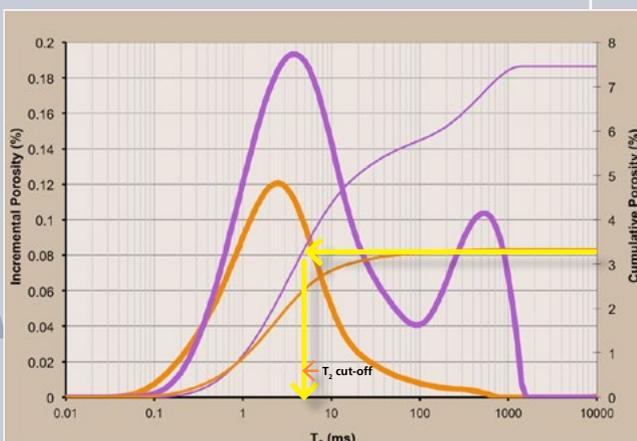


Figure 2 - T_2 distributions showing T_2 cut-off point

An accurate determination of the T_2 cut-off point is essential for an accurate determination of recoverable reserves (mobile fluid). T_2 cut-off can be easily determined in the laboratory by using two NMR measurements; one on a cleaned and re-saturated plug, the other on the same plug after it has been spun in a centrifuge to irreducible water saturation. T_2 distributions are plotted for both data sets, along with the cumulative values of the distributions. The T_2 cut-off is taken to be the point at which the cumulative value of the saturated distribution (yellow horizontal arrow in Figure 2) equals the final cumulative value of the irreducible distribution (vertical yellow arrow in Figure 2). The data plotting and calculation for this measurement are carried out automatically by the LithoMetrix software supplied as standard with every **GeoSpec** NMR core analyser.

This measurement should be done on core plugs taken from various zones in the well, with the resulting T_2 cut-off values applied retrospectively to the log data from each depth. This will result in a more accurate prediction of Bound Volume Index (BVI) and Free Fluid Index (FFI) than those obtained by applying a single default value of the T_2 cut-off to the entire log.

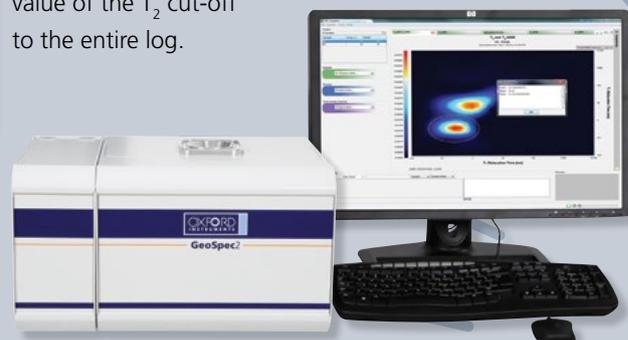

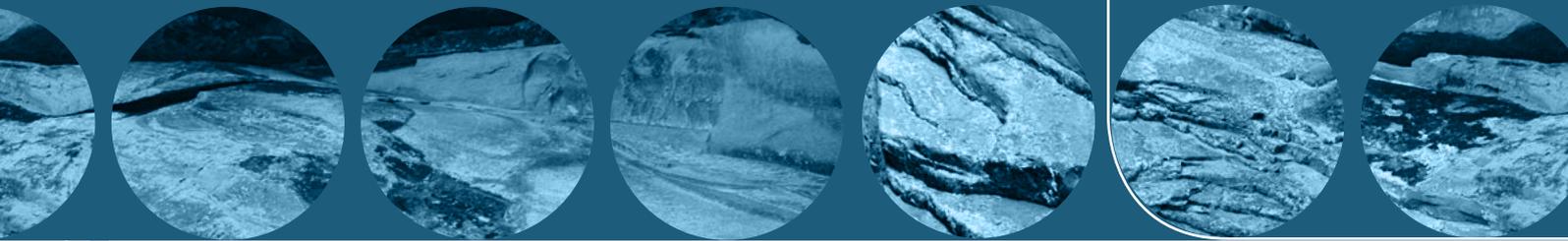

The following example illustrates how reserve estimates can be seriously under-reported if a standard cut-off value is used without laboratory calibration. The samples used in this study were from a tight sandstone play.

In this particular example, laboratory measurements correct the standard cut-off estimate of free fluid index by an average of 260%. This means that this reserve could produce 2.6 times more oil than would have been estimated if the laboratory calibration had not been completed.

Sample	Porosity (%)	Permeability (mD)	FFI at 33msec (%)	Perm. Est. (mD)	Measured FFI (%)	Cut-off (ms)
1	5.71	0.003	0.46	0.000068	1.34	10.00
2	7.59	0.003	1.04	0.000080	2.78	7.08
3	4.59	0.001	0.24	0.000014	0.47	14.13
4	7.36	0.003	1.46	0.00170	3.13	7.94
5	12.01	0.28	2.71	0.045	7.88	1.79
6	10.03	0.22	1.51	0.004	2.42	10.66
7	10.2	0.13	2.00	0.032	4.42	7.77
8	8.8	0.15	1.17	0.031	5.08	2.66

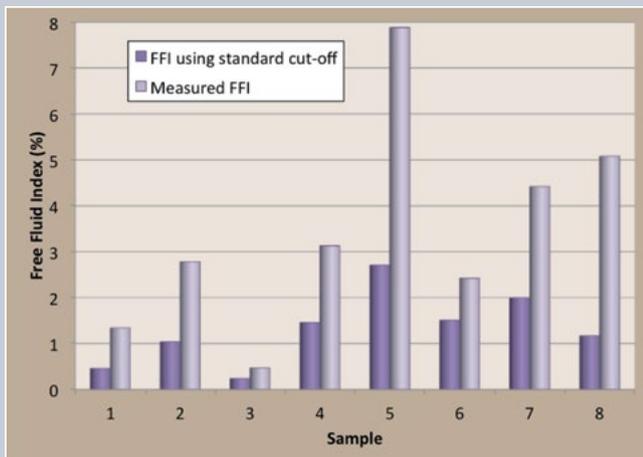


Figure 3 – Free Fluid Index using standard and measured cut-off values

Permeability

Permeability can be calculated from the T_2 distribution data using one of two commonly accepted mathematical models: the free-fluid or Coates model (Equation 1) can be applied in formations containing water and/or hydrocarbons, while the average- T_2 or Schlumberger model (Equation 2) can be applied to pore systems containing only water and gas. In either case, measurements on core samples are necessary to refine these models by determining the correct values of the coefficients, and produce a model customised for local use.

$$\kappa = \left(\frac{\phi}{c}\right)^a \left(\frac{\text{FFI}}{\text{BVI}}\right)^b \quad \kappa = aT_2^2 \log \text{mean} \phi^4$$

Equation 1 - The Coates model for permeability

Equation 2 - The Schlumberger model for permeability

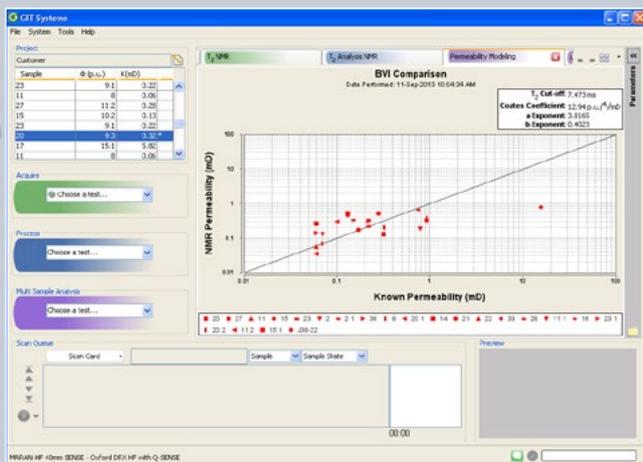


Figure 4 –NMR permeability versus actual

Customisation of these models is done by inserting values of the log mean T_2 , FFI and BVI obtained from core plug NMR measurements, and then varying the coefficients until the best linear fit is obtained against permeability measurements determined by other techniques. **LithoMetrix** (included with every **GeoSpec**) provides tools to perform this on single samples while **LithoMetrix** has advanced regression tools allowing optimum T_2 cut-offs and coefficients to be determined for a set of samples. Figure 4 shows NMR permeability versus actual permeability using the optimized values for a large set of samples.

GeoSpec – The new approach to core analysis using NMR



Necessary apparatus

Tools needed to perform NMR log calibrations as described in this note include: a **GeoSpec** NMR Core Analyser with **LithoMetrix** software; a centrifuge capable of removing all mobile (free) fluid from the core plug; and equipment to measure permeability. Although not required to perform log calibrations, pulsed field gradients on the **GeoSpec** instrument and optional **GIT Systems** software can allow users to perform additional advanced measurements to further enhance interpretation of NMR logs.

Recommendations for consideration

Expanding our supply of oil and gas requires exploration. With increasing focus on unconventional reserves, it is imperative that we do not rely on previously acceptable “standards” without question. Example after example shows that using standard cut-offs significantly misrepresents reserve estimates. Easy-to-use **GeoSpec** NMR systems allow rapid, consistent, reliable laboratory measurements to support work completed at the well site; providing reliable results for the reservoir engineers and petrophysicists.



Oxford Instruments Industrial Analysis

For more information please email:

industrial@oxinst.com

www.oxford-instruments.com

UK

Tubney Woods, Abingdon, Oxfordshire, OX13 5QX, UK

Tel: +44 (0) 1865 393 200 **Fax:** +44 (0) 1865 393 333

USA

300 Baker Avenue, Suite 150, Concord, MA, 01742, USA

Tel: +1 978 369 9933 **Fax:** +1 978 369 8287

China

Floor 1, Building 60, No. 461 HongCao Road,
Shanghai, 200233, PRC, China

Tel: +86 21 61273858 **Fax:** +86 21 6360 8535

Green Imaging Technologies

For more information: info@greenimaging.com

www.greenimaging.com

Canada

520 Brookside Drive, Suite B, Fredericton, NB, E3A 8V2,
Canada

Toll Free: +1 888 944 8462

Tel: +1 506 458 9992 **Fax:** +1 506 458 9615

This publication is the copyright of Oxford Instruments plc and provides outline information only, which (unless agreed by the company in writing) may not be used, applied or reproduced for any purpose or form part of any order or contract or regarded as the representation relating to the products or services concerned. Oxford Instruments' policy is one of continued improvement. The company reserves the right to alter, without notice the specification, design or conditions of supply of any product or service. Oxford Instruments acknowledges all trademarks and registrations. © Oxford Instruments plc, 2015. All rights reserved. Ref: LogCal-AN-07-15



The Business of Science®